


LONGCIP

LONGCIP-1, LONGCIP-2-3 Product Description

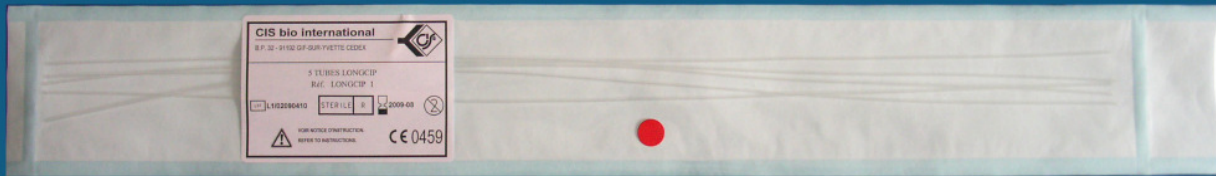

BEBIG LONGCIP are used for interstitial, flexible implants in temporary brachytherapy, e.g. head and neck, breast and soft tissue implants, etc.


■ LONGCIP-1

LONGCIP-1 is a flexible implant tube and serves as a carrier for LONGCIP-2-3. LONGCIP-2-3 is inserted into LONGCIP-1 by using a leader insertion needle.


■ LONGCIP-2-3

LONGCIP-2-3 is designed to guide BEBIG Ir-192 Wire Sources which have a diameter of 0.3 mm. The Ir-192 Wire Source is directly inserted into LONGCIP-2-3. LONGCIP-2-3 can also be used as leader for LONGCIP-1.


■ Treatment Examples

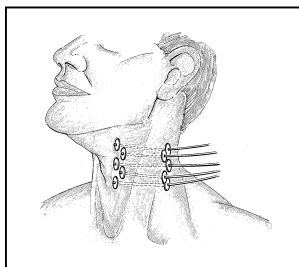

Neck

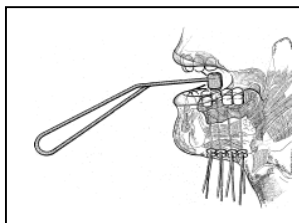

Tongue

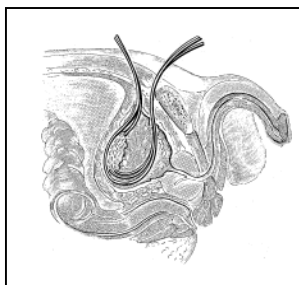


Bladder

Distributed by SeedOS Ltd
Please contact cwalters@seedos.com
www.seedos.com

 Eckert & Ziegler
BEBIG

LONGCIP


■ LONGCIP Details

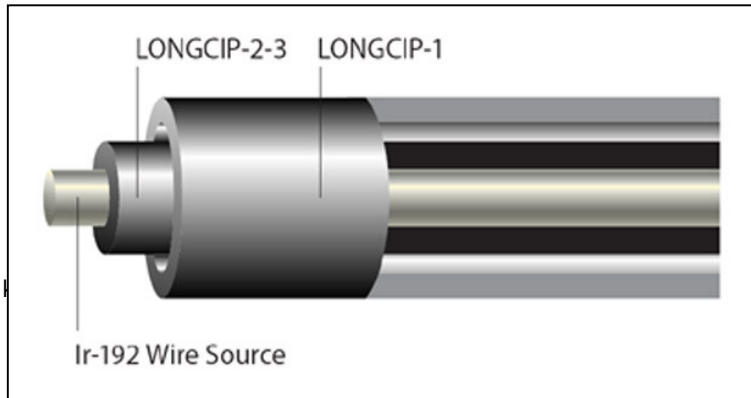

■ Technical Data

	LONGCIP-1	LONGCIP-2-3
Length	500 mm	500 mm
Inner Diameter	1.2 mm	0.45 mm
Outer Diameter	1.6 mm	0.90 mm

■ Flexible implant loop insertion technique for the base of the tongue

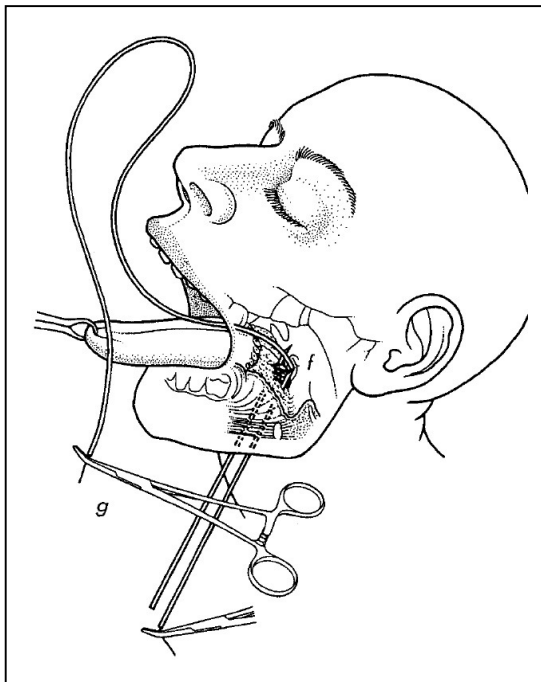

Hilaris, Nori, Anderson, An Atlas of Brachytherapy, New York, 1988, p. 127

**Eckert & Ziegler
BEBIG GmbH**

Robert-Rössle-Strasse 10
13125 Berlin
Germany
www.bebig.com

Telephone +49 (0) 30 94 10 84-130
Telefax +49 (0) 30 94 10 84-112
e-mail bebig@ezag.de